
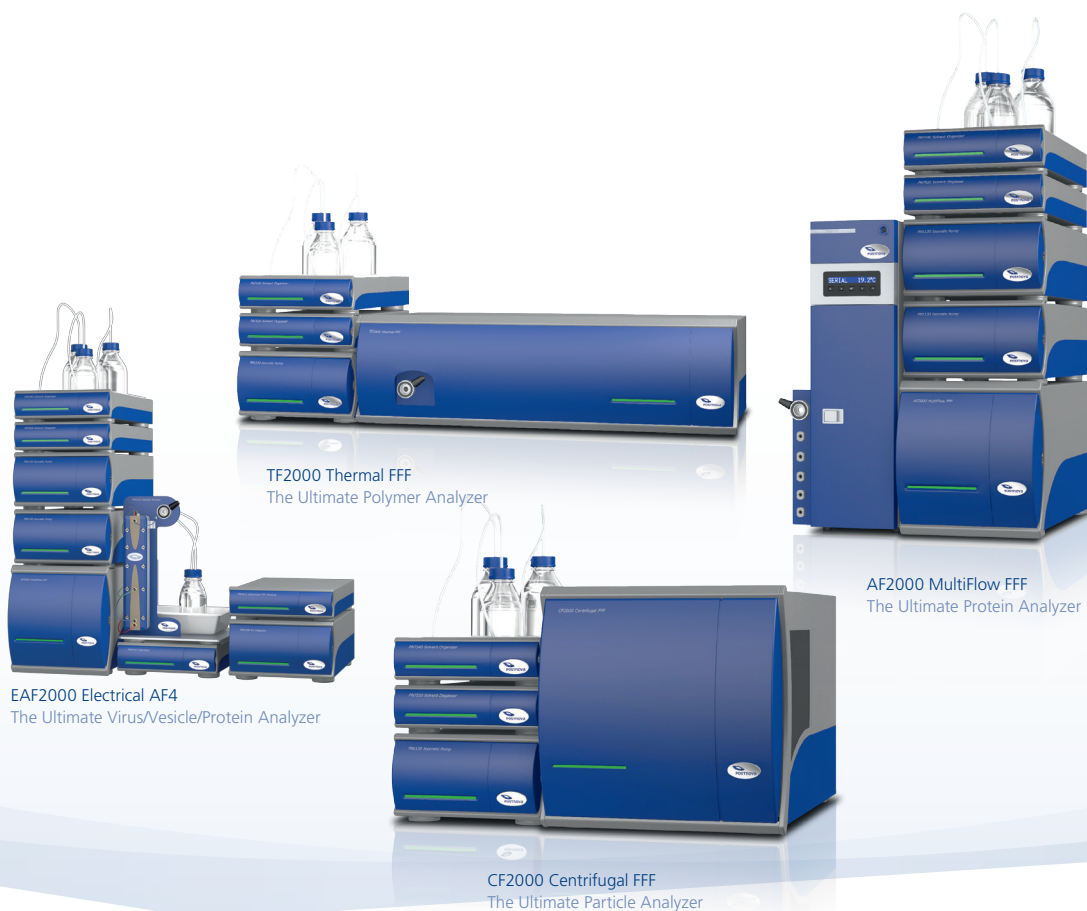


NovaFFF Software Platform

NovaFFF Control Software Platform for Postnova FFF Systems

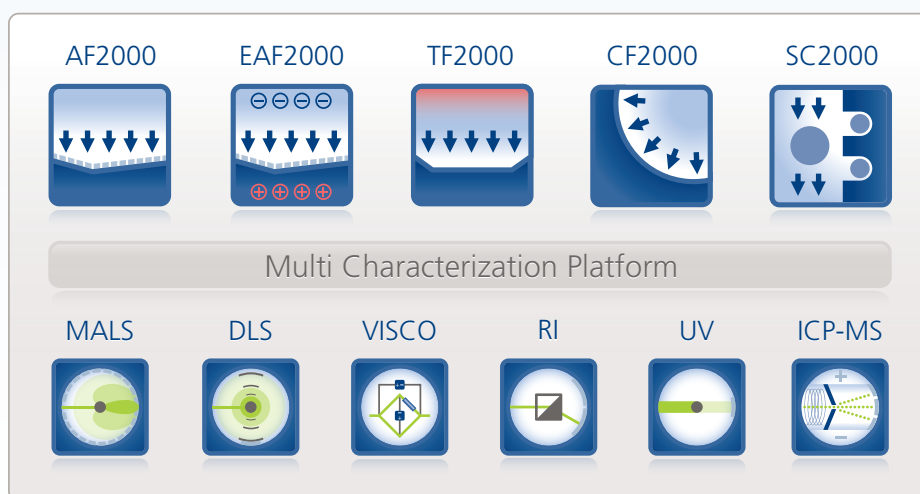

NovaFFF Software Platform

NovaFFF is an advanced software platform suite for Postnova's Field-Flow Fractionation instruments. It controls instruments, collects data, and allows users to evaluate data and generate reports.

Features

- **Analysis and Calculation of FFF Results**

Import and analysis of FFF data

Data processing of data from Asymmetrical Flow FFF (AF4), Electrical AF4 (EAF4), Centrifugal FFF (CF3) and Thermal FFF (TF3)

Calculation of particle size, distribution and molar mass based on published FFF theory

- **Easy to Use - the NovaFFF Software**

Controls the FFF platform together with the MALS detector

Automatic system shutdown if errors detected (e.g. over-pressure)

Simple integration of additional hardware components e.g. autosampler, fraction collector and extra detectors

Continuously monitors the system pressures, including pumps and channel


Delivers comprehensive raw data files that include pressures, flows and all detector signals

Outputs a log file providing full data file history

Ordering Information

P-NVA-FFF-001

NovaFFF Software

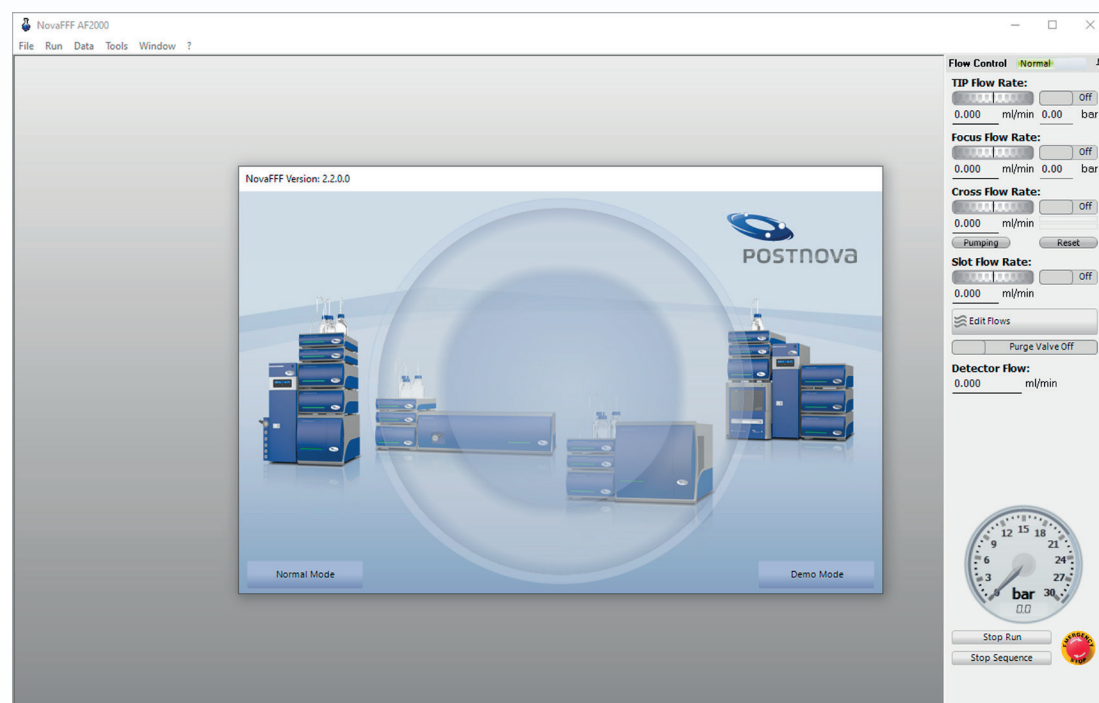

Specifications

- Supported Operating Systems:
Windows 10 – 64 bit
Windows 11 – 64 bit (*)
(*recommended)
- Hardware Requirements:
Intel® Core® i5 – 7500 (3.41 GHz)
recommended Intel® Core® i7 – 7500 or better

16 GB RAM or more

Color Display min. 1920 x 1080

Contact

- Postnova Analytics GmbH
86899 Landsberg, GERMANY
T: +49 8191 985 688 0
- Postnova Analytics UK Ltd.
Malvern, Worcestershire, WR14 3SZ, UK
T: +44 1684 585167
- Postnova Analytics Inc.
Salt Lake City, UT 84102, USA
T: +1 801 521 2004

info@postnova.com
www.postnova.com